



LVS Series

Large Vet Scales



LVS 700 / LVS XL 700



Reliable!

LVS + 700 / LVS XL + 700



Accurate!



LVS-S 700

Quality

LW Measurements LLC takes pride in producing quality products that are affordable, durable, and most importantly accurate. Visit our website for more information and videos of our products.

www.lwmeasurements.com



Large Vet Scale

The Vet and Agriculture community love our vet scale series. This series is ideal for weighing dogs, exotic animals, livestock, goats, pigs and more. The large stainless steel weighing platform was crafted to withstand daily use. Featuring an extended digital read out with two modes of power, AC Adapter and Battery options.

Features

The LVS series is manufactured with top-of-the-line software and hardware offering the best product without the hefty price tag. The three variations of the LVS 700 models provided have a max capacity of 700lb x 0.2lb. Every LVS has a steel frame and extra-large stainless-steel base/ platter, with structural integrity for the larger animals and easy sanitation. These scales are intentional in their design to help provide the most efficient and accurate weighing.

LVS 700 and XL 700

Choosing the correct LVS made simple. Both offer a 700lb max capacity, stainless steel platter, wired REMO indicator, adjustable feet, caster wheels, a handle, and a rubber protective mat. What separates the LVS and the LVS XL is the platter size. The LVS XL offers 6 more inches in length and 3 inches in width.



LVS + Model

The LVS + Model offers the same capacity, platter size, and features as the LVS 700. What separates this from its counterpart is the Bluetooth capability. The Bluetooth feature offers animal safety, storage, and portability. The scale was designed after receiving calls of multiple animals chewing the cords. Every (+) model has a Bluetooth indicator and Rechargeable lithium-ion batteries in the scale base and indicator. Removing the stress of being tethered to an electrical socket.



LVS-S 700. This model is the stowaway/compact version of the LVS 700 products. It has different profile since it is mounted the wall. Which allows you to fold it up and out of the way after every use. Features dual power source: AC power adapter or batteries.



Convenience

Each scale includes AC power adapters, rubber protective mats, adjustable feet, and wheels (LVS regular and plus models only). The LVS 700 comes with a Remo indicator that is factory set for the specific capacity of this scale model. The LVS 700 models have two modes of power allowing for portable use. A RS232C cable port is included to connect to computers and printers. The LVS-S is our wall mounted option of this model, which allow for easy stowing when not in use. The stowaway is a perfect option for vet clinics that have limited space.

The Remo Indicator



User friendly interface.

Hold, Zero, Tare, and Print at the press of a button. With the ability to connect a printer or to a computer.

Battery Compartment



AC Adapter Port

RS232C



Versatility

A durable handle and wheels have been added to make transferring or moving the scale easier. Adjustable footing was designed for areas that have uneven surfaces to create a safe and stable platform for the animal being weighed.






All LVS models come with a Rubber Protective Mat



LVS-S wall mounted space saver ...

The LVS-S was designed to be mounted on a wall and stored up and out of the way after use to help save space and convenience.

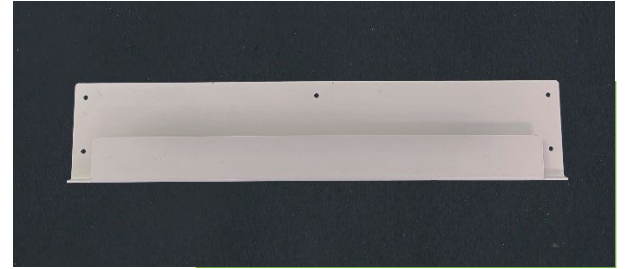

Features

LVS 700 & LVS XL 700

- 700lb x 0.2lb capacity
- REMO indicator
- Auto shut off
- AnyCal Calibration
- Auto backlight function
- Units: lbs. / kg (option)
- Stainless platter
- LVS 700 platter size 37.4" x 19.7"
- LVS XL 700 platter **size** 44" x 22"
- AC adapter or 6 x AA batteries

LVS + 700 & LVS XL + 700 (wireless)

- 700lb x 0.2lb capacity
- Wireless REMO indicator and scale base
- Auto shut off
- AnyCal Calibration
- Auto backlight function
- Units: lbs. / kg (option)
- Stainless platter
- LVS +700 platter size 37.4" x 19.7"
- LVS XL +700 platter **size** 44" x 22"
- AC adapter or Rechargeable lithium-ion batteries

LVS-S 700

- 700lb x 0.2lb capacity
- REMO indicator
- Wall Mounted stowaway
- Auto shut off
- AnyCal Calibration
- Auto backlight function
- Units: lbs. / kg (option)
- Stainless platter
- LVS-S 700 platter size 37.4" x 19.7"
- AC adapter or 6 x AA batteries

Specifications



	LVS 700	LVS XL 700	LVS-S 700
Model Suffix	LVS 700 / LVS + 700	LVS XL 700 / LVS XL + 700	LVS-S 700
Casing	Steel		
Capacity	700lb		
Resolution(kg)	0.1		
Resolution(lb.)	0.2		
Resolution(oz)	No		
Display Divisions	3500		
Internal Resolutions	2 m		
Accuracy	+/-3d		
Platter Size(inch)	37.4 x 19.7	44 x 22	37.4X19.7
Scale Dimensions(inch)	37.4 x 19.7 x 3	44 x 22 x 3	37.4X19.7X2.8
LCD Height(inch)	1		
Beeper	Yes		
Auto Backlight	Yes, Selectable		
Auto Shut Off	Yes, Selectable		
Auto Hold	No		
Unit Switching	kg, lb.		
AD Processor	21-bit A/D		
AD Update Speed	10 Hz		
Settling Time	<2s		
Calibration Software	AnyCAL		
No of Keys	5		
Memory	Zero		
Counting Function	No		
RS232C	Yes		
USB	No		
Battery	Standard and stowaway models use 6 x AA Batteries / All (+) Models use rechargeable lithium-ion batteries		
AC Adapter	9V 400mah		
Operating Environment	0°C - 40°C / 32°F - 104°F		
Standard Accessories	Manual, AC adapter, Rubber protective mat		
Accessories sold separately	All LVS models have replaceable indicators		

WARRANTY POLICY

LW Measurements products are under warranty against factory defects for a period of two (2) years from the date of shipment for Customers in the continental United States.

LW Measurements will pay for freight both ways within the first 30 days of purchase if technical support is required. (Determined by our certified Technical Support Team) In the event technical support is needed the customer is responsible for shipping the product back to us after the 30-day period has expired. Once the product is received, the technical support team will inspect and as necessary repair or replace the scale. LW Measurements will ship the product back to the customer at our expense.

Any new scales returned for warranty must be properly packaged in the original box. If they are not properly packed in the original box, the customer is responsible for the shipping cost. If we determine there is a factory defect, we will provide a replacement and accept responsibility for shipping back to customer. If we determine that it is not a factory defect, the customer will be responsible for shipping fees incurred.

For Customers outside the lower 48 States, including Mexico, Canada, Puerto Rico, Hawaii, Alaska and all other countries are responsible for all shipping fees incurred.

Our warranty does not cover misuse or neglect including but not limited to battery or water damage, overloading, and chewed or cut wires, keypad damage (punctures), scratched platforms, dented platforms unless caused by shipping (Must report/return within 30 days). If the product is found to have been misused or damaged by the customer, LW measurements is not responsible for the cost of return.

For warranty claims, please call your customer service representative.

